

TS 031 – CHROMOGENIC ENRICHMENT SUPPLEMENT

INTENDED USE

For the isolation, detection and enumeration of *Listeria monocytogenes*.

COMPOSITION

Ingredients	Concentration
L – phosphatidylinositol	10ml

(Per vial sufficient for 500 ml medium)

Applied in Medium

TM 004 - CHROMOGENIC LISTERIA AGAR BASE (ISO 11290-1 & 2:2004, ISO 11133:2014)


TM 1845 - CHROMOGENIC L-MONO LISTERIA DIFFERENTIAL AGAR

INSTRUCTION FOR USE

Aseptically add the content of 2 vials of CHROMOGENIC ENRICHMENT SUPPLEMENT (TS 031) at slightly high temperature (70°C) with vigorous shaking to 970ml of sterile molten CHROMOGENIC LISTERIA AGAR BASE (TM 004) / CHROMOGENIC L-MONO LISTERIA DIFFERENTIAL AGAR (TM 1845) along with sterile rehydrated contents of 2 vial of CHROMOGENIC LISTERIA SELECTIVE SUPPLEMENT (TS 205). Mix well and pour into sterile Petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

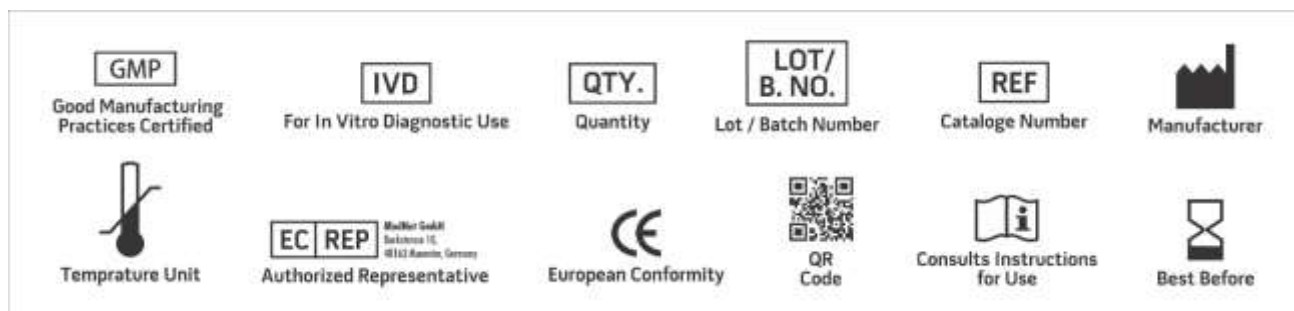

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 29 JULY 2024