

1205V-VEG. B. MEAT EXTRACT POWDER

INTENDED USE

B. Meat extract paste (Bovine) are used in preparing microbiological culture media

PRODUCT SUMMARY AND EXPLANATION

B. MEAT EXTRACT POWDER (Type-1) is intended to replace aqueous infusion of meat in microbiological culture media. It is frequently used at a concentration of 0.3 to 1.0% in culture media, although concentrations may vary depending on the nutritional requirements for the medium formulation. It may be relied upon for biochemical studies, particularly fermentation reactions, because of its independence from fermentable substances that would interfere with the accuracy of such determinations.

PRINCIPLE

B. MEAT EXTRACT POWDER (Type-1) is infusion of beef and provide an undefined source of nutrients. Beef Extract products are not exposed to the harsh treatment used for protein hydrolysis, so they can provide some of the nutrients lost during peptone manufacture. Beef Extract Powder is mixtures of peptides and amino acids, nucleotide fractions, organic acids, minerals and some vitamins. The function of B. Meat extract paste (Bovine) products can be described as complementing the nutritive properties of peptone by contributing minerals, phosphates, energy sources and those essential factors missing from peptone.

INSTRUCTION FOR USE

B. MEAT EXTRACT POWDER (Type-1) is a dehydrated extract of bovine tissue for use in preparing microbiological culture media in a laboratory setting. Beef Extract Powder is not intended for use in the diagnosis of disease or other conditions in humans. It is prepared and standardized for use in microbiological culture media, where it is generally used to replace infusion of meat. Culture media containing B. Meat extract powder Powder are recommended for use in bacteriological examination of water, milk, and other materials, where uniform composition of media is important. It is relied upon for biochemical studies, particularly fermentation reactions because of its independence from fermentable substances.

QUALITY CONTROL SPECIFICATIONS

Appearance	:	Light yellowish to brownish yellow color, free flowing
Solubility (2% soln. at 25°C)	:	Soluble in distilled water, clear.
Clarity (2% Soln. at 121°C)	:	Clear solution. No ppt.
pH (2% Soln. at 25°C)	:	6.5 – 7.5
Loss on drying (at 105°C)	:	NMT – 6.0%
Total Nitrogen (DWB)	:	NLT – 12.0%
α-Amino Nitrogen	:	NLT – 2.5%
Total Ash	:	NMT – 12.0%
Chloride (as NaCl)	:	NMT – 6.0%
Indole Test	:	Positive
Microbial Test	:	Passes Test

INTERPRETATION

Cultural Characteristic observed in 2% Veg. BE Extract Powder and 1.5% agar after incubation at 35-37°C for 18-24 hours

Microorganism	ATCC	Inoculum (CFU/ml)	Growth
<i>Staphylococcus aureus</i>	25923	50-100	Good - Luxuriant
<i>Escherichia coli</i>	25922	50-100	Good - Luxuriant



<i>Pseudomonas aeruginosa</i>	27853	50-100	Good - Luxuriant
<i>Bacillus subtilis</i>	6633	25923	Good - Luxuriant
<i>Enterococcus faecalis</i>	29212	25922	Good - Luxuriant
<i>Streptococcus pyogenes</i>	19615	27853	Good - Luxuriant

PACKAGING

Standard packing is 500gm in HDPE bottle. After packing tightly closed in a dry and well-ventilated place.

STORAGE

Store at room temperature in cool place, Keep container tightly closed in a dry and well-ventilated place and away from bright light. Use before expiry date on label. On opening, product should be properly stored in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use.

Product Deterioration: Do not use product if any contamination, discoloration or other sign of deterioration is found.

DISPOSAL


After use, contact a licensed professional waste disposal service to dispose of this material. Dispose of as unused product.

REFERENCES

1. Prokofeva, Miroshnichenko, Kostrikina, Chernyh, Kuznetsov, Tourova and Bonch-Osmolovskaya. 2000. Int. J. Syst. Evol. Microbiol. 50: Pt 6:2001.
2. United States Pharmacopeial Convention, Inc. 2001. The United States pharmacopeia 25/The national formulary 20 – 2002. United States Pharmacopeial Convention, Inc., Rockville, Md.
3. Clesceri, Greenberg and Eaton (ed.). 1998. Standard methods for the examination of water and wastewater, 20th ed. American Public Health Association, Washington, DC.
4. U.S. Food and Drug Administration. 1995. Bacteriological analytical manual, 8th ed. AOAC International, Gaithersburg, Md.
5. Downes and Ito (ed.). 2001. Compendium of methods for the microbiological examination of foods, 4th ed. American Public Health Association, Washington, D.C.
6. Horwitz (ed.). 2000. Official methods of analysis of AOAC International, 17th ed. AOAC International, Gaithersburg, Md.
7. U.S. Department of Agriculture. 1998. Microbiology laboratory guidebook, 3rd ed. Food Safety and Inspection Service, USDA, Washington, D.C.
8. Cote. 1999. In Flickinger and Drew (ed.), Encyclopedia of bioprocess technology: fermentation, biocatalysis, and bioseparation. John Wiley & Sons, Inc., New York, N.Y. 10. Bridson and Brecker. 1970. In Norris and Ribbons (ed.), Methods in microbiology, vol. 3A. Academic Press, New York. N.Y.



Quantity

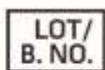

Lot / Batch Number



Temperature Unit



Best Before



QR Code



Catalogue No.



Consults Instructions for use :



Manufacturer

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**

Revision: 05th Oct. 2019

