

1581 – PEPTONE TYPE-II (Bacteriological Grade)

INTENDED USE

Peptone Type-II used in the preparation of culture media employed for cultivation of a wide variety of microorganisms.

PRODUCT SUMMARY AND EXPLANATION

Peptone Type-II is used in preparing microbiological culture media and in producing bacterial toxins and also usable in synthetic media in acclimatization of microorganisms in bioreactor studies. It's support to growth of Staphylococci, Streptococci, Pneumococci and also suitable for isolating and cultivating Haemophilus and Neisseria. It is off white to Creamish yellow colour, free flowing powder having characteristic odour but not pungent smell. It is completely soluble in distilled Water, Clear. Insoluble in alcohol.

PRINCIPLE

Peptone Type-II is enzymatic digest of protein used in preparing microbiological culture media and in producing bacterial toxins. Proteose peptone provide nitrogen in a form that is readily available for bacterial growth. It is superior in nutritious of fastidious microorganism.

INSTRUCTION FOR USE

Peptone Type-II is used in media for the production of bacterial toxins. It is used in preparing chocolate agar for propagating of Neisseria species. It is also used for the cultivation of bacteria with high nutritional requirements, as for example Haemophilus, Salmonella, staphylococcus etc. species.

QUALITY CONTROL SPECIFICATIONS

Appearance	:	Light yellowish to brownish yellow colour, free flowing powder having characteristic odour but not pungent smell.
Solubility (2% Soln. at 25°C)	:	Completely soluble in distilled Water, Clear. Insoluble in alcohol.
pH (2% Soln. at 25 °C)	:	6.5 – 7.5
Loss on drying (at 105 °C)	:	NMT – 7.0%
Total Nitrogen (DWB)	:	NLT – 13.5%
α-Amino Nitrogen	:	NLT – 3.5%
Total Ash	:	NMT – 12.0%
Chloride (as NaCl)	:	NMT – 5.0%
Indole Test	:	Positive
Microbial Parameter	:	Passes Test
Growth Promotion Test	:	Passes Test

INTERPRETATION

Cultural Characteristic observed in 2% Peptone Type-III and 1.5% agar after incubation at 35-37°C for 18-24 hours.

Microorganism	ATCC	Inoculum (CFU/ml)	Growth
Staphylococcus aureus	6538	50-100	Luxuriant
Escherichia coli	8739	50-100	Luxuriant
Pseudomonas aeruginosa	9027	50-100	Luxuriant
Salmonella typhi	6539	50-100	Luxuriant
Streptococcus pyogenes	19615	50-100	Luxuriant



PACKAGING:

Standard packing is 500gm in plastic bottle. After packing tightly closed in a dry and well-ventilated place.

STORAGE

Keep plastic bottle tightly closed in a dry and well-ventilated place, Store in cool place. Use before expiry date on label. On opening, product should be properly stored in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the plastic bottle after use.

Product Deterioration: Do not use product if any contamination, discoloration or other sign of deterioration is found.

DISPOSAL


After use, contact a licensed professional waste disposal service to dispose of this material. Dispose of as unused product.

REFERENCES

1. Kirkbride, Berthelsen and Clark. 1931. Comparative studies of infusion and infusion-free diphtheria toxin in antitoxin production and in standardization by the flocculation, subcutaneous, and intracutaneous tests. J. Immunol. 21:1-20.
2. Hazen and Heller. 1931. Further studies upon the effect of various carbohydrates on production of diphtheria toxin with special reference to its flocculating titer and final pH. J. Bacteriol. 23:195-209.



Quantity

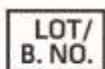

Lot / Batch Number



Temperature Unit



Best Before



QR
Code



Catalogue No.



Consults Instructions for use :



Manufacturer

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 05th Oct. 2019