

TS 007 – CAMPYLOBACTER SUPPLEMENT-I (Blaser Wang)

INTENDED USE

For selective isolation of *Campylobacter* species.

COMPOSITION

Ingredients	Concentration
Polymyxin B	1.250IU
Vancomycin	5mg
Trimethoprim	2.500mg
Amphotericin B	1mg
Cephalothin	7.500mg

(per vial sufficient for 500 ml medium)

Applied in Medium

TM 041 - BLOOD AGAR BASE No. 2
 TM 1154 - C. BOTULINUM ISOLATION AGAR BASE
 TM 056 - CAMPYLOBACTER AGAR BASE
 TM 071 - COLUMBIA BLOOD AGAR BASE
 TMV 071 - COLUMBIA BLOOD AGAR BASE (VEG.)
 TM 1842 - BLOOD AGAR BASE NO. 2 (ISO 11290-1:2017)
 TM 2011 - BLOOD AGAR BASE NO. 2 WITH 1.2% AGAR
 TM 514 - BRUCELLA AGAR BASE
 TMV 514 - BRUCELLA AGAR BASE (VEG.)
 TM 050 - BRUCELLA BROTH BASE
 TM 1157 - CAMPYLO THIOGLYCOLLATE MEDIUM BASE (CAMPY-THIO MEDIUM)
 TM 2044 - COLUMBIA BLOOD AGAR BASE W/ 1% AGAR


INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml sterile distilled water and mix well to dissolve. Avoid frothing of the solution. Aseptically add the contents to 500 ml of sterile, molten, cooled (45-50°C) TM 041 - Blood agar base no. 2 / TM 1154 - C. Botulinum isolation agar base / TM 056 - Campylobacter agar base / TM 071 - Columbia blood agar base / TMV 071 - Columbia blood agar base (veg.) / TM 1842 - Blood agar base no. 2 (ISO 11290-1:2017) / TM 2011 - Blood agar base no. 2 with 1.2% agar / TM 514 - Brucella agar base / TMV 514 - Brucella agar base (veg.) / TM 050 - Brucella broth base / TM 1157 - Campylo thioglycollate medium base (campy-thio medium) / TM 2044 - Columbia blood agar base w/ 1% agar with 10% v/v sheep blood or 5-7% v/v laked horse blood and rehydrated contents of 1 vial of Campylobacter Growth Supplement. Mix gently and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.



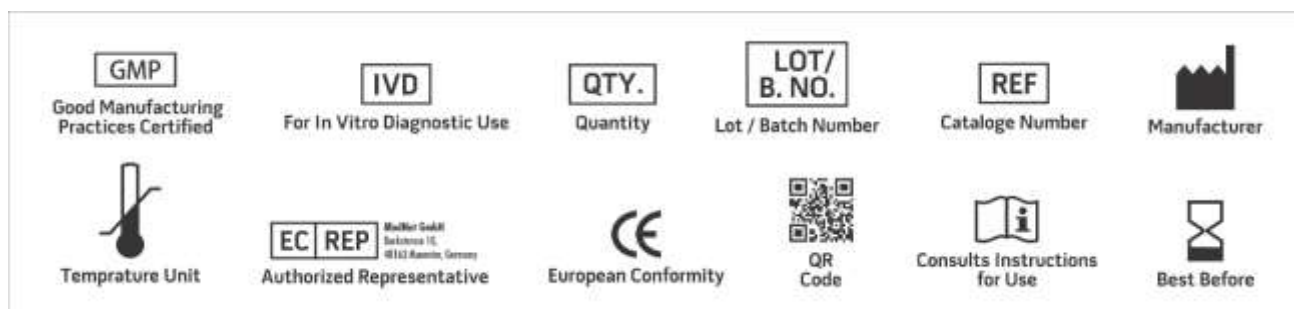

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 11 Feb., 2022