

## TS 028 – DERMATO SUPPLEMENT

### INTENDED USE

For the determine the tryptophan-deaminase activity.

### COMPOSITION

Ingredients	Concentration
Cycloheximide	250mg
Chlortetracycline	50mg
Gentamicin	50mg

(per vial sufficient for 500 ml medium)

### Applied in Medium


TM 083 - D.T.M. AGAR BASE (DERMATOPHYTE TEST AGAR BASE)

### INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 5 ml of 50% v/v acetone. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) TM 083 - D.T.M. AGAR BASE. Mix well and pour into sterile petri plates.

### STORAGE

Vials should be stored in sealed container at 2-8°C.

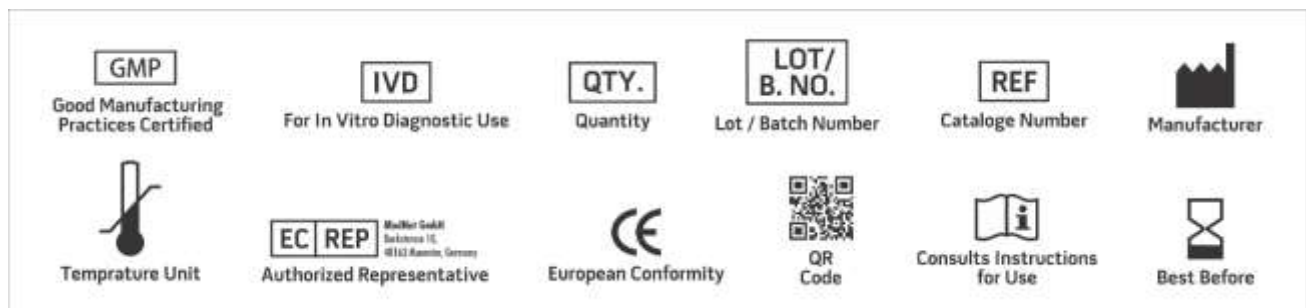

**NOTE:** Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

**\*For Lab Use Only**  
Revision: 28 Feb., 2022