

TS 045 – BROMO CRESOL PURPLE (15 mg/vl)

INTENDED USE

For selective isolation and enumeration of faecal Streptococci.

COMPOSITION

Ingredients	Concentration
Bromo cresol purple	15mg

(Per vial sufficient for 1000 ml medium)

Applied in Medium


TM 136 - KF STREPTOCOCCAL AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml of sterile 0.2 N sodium hydroxide and add to 1000 ml of sterile, molten TM 136 - KF STREPTOCOCCAL AGAR BASE, along with 1 vial of 1% TTC Solution. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

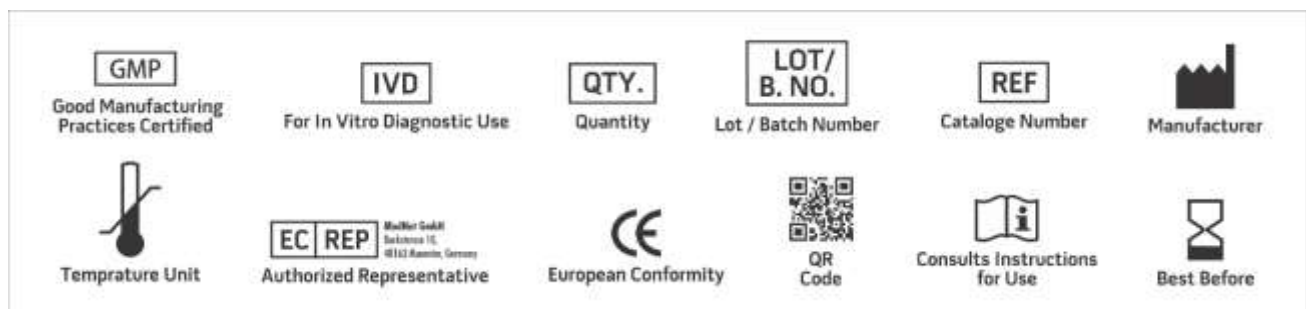

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 5 March., 2022