

TS 068 – KANAMYCIN SULPHATE SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of Enterococci.

COMPOSITION

Ingredients	Concentration
Kanamycin sulphate	10mg

(per vial sufficient for 500 ml medium)

Applied in Medium

TM 748 - KANAMYCIN ESCULIN AZIDE AGAR BASE


TM 749 - KANAMYCIN ESCULIN AZIDE BROTH BASE

INSTRUCTION FOR USE

Rehydrate the content of one vial aseptically with 2 ml of sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) TM 748 - KANAMYCIN ESCULIN AZIDE AGAR BASE / TM 749 - KANAMYCIN ESCULIN AZIDE BROTH BASE. Mix well and aseptically pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

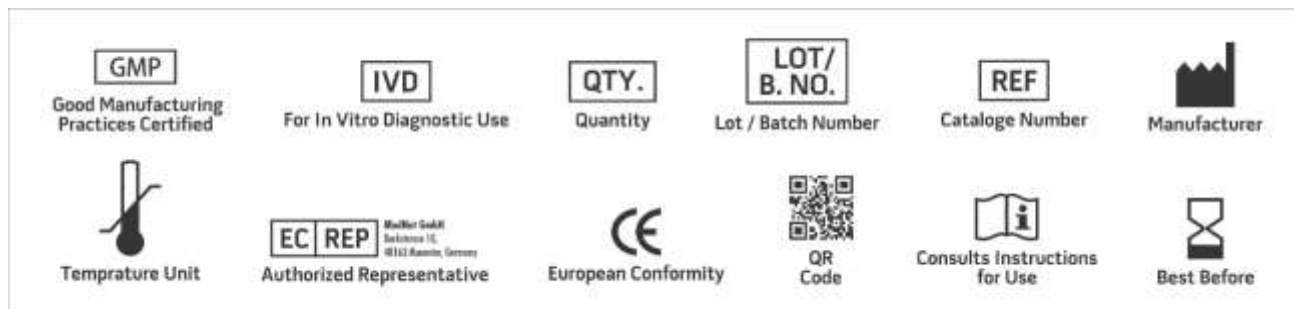

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 2 March., 2022