

TS 072 – XLT4 SUPPLEMENT

INTENDED USE

For isolation of non-typhi *Salmonella*.

COMPOSITION

Ingredients	Concentration
(40ml per vial contains Tergitol 4 - 27%)	-
Tergitol 4	10.800g
Distilled water	40ml

(per vial sufficient for 1000 ml medium)

Applied in Medium


TM 493 - XLT4 AGAR BASE

INSTRUCTION FOR USE

Warm up XLT4 Supplement to room temperature and add 4.6ml to 1000 ml of TM 493 - XLT4 AGAR BASE. Heat to boiling to dissolve the medium completely. Avoid overheating. DO NOT AUTOCLAVE. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

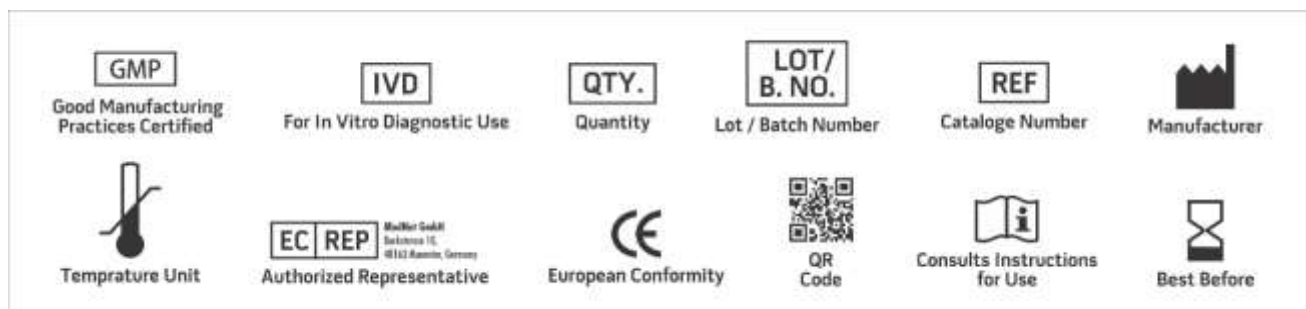

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only
Revision: 2 March., 2022