

TS 075 – CETRINIX SUPPLEMENT

INTENDED USE

For the selective isolation of *Pseudomonas* species.

COMPOSITION

Ingredients	Concentration
Cetrimide	100mg
Nalidixic acid	7.500mg

(Per vial sufficient for 500 ml medium)

Applied in Medium


TM 626 - PSEUDOMONAS AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2 ml of sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) TM 626 - PSEUDOMONAS AGAR BASE. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

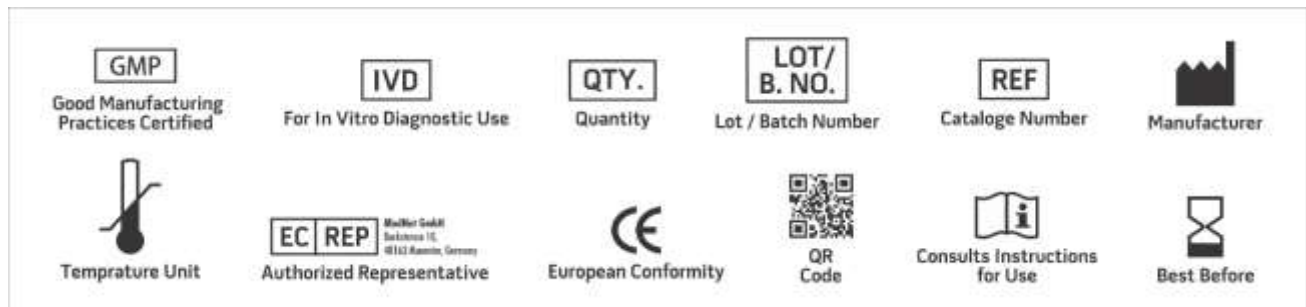

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 8 March., 2022