

TS 100 – BACTEROIDES SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of *Bacteriodes* species.

COMPOSITION

Ingredients	Concentration
Gentamicin	50mg

(Per vial sufficient for 500 ml medium)

Applied in Medium


TM 1144 - BACTEROIDES BILE ESCULIN AGAR BASE (BBE AGAR BASE)

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 5.0 ml sterile distilled water. Mix well and aseptically add to 500 ml sterile molten, cooled (45-50°C) TM 1144 - BACTEROIDES BILE ESCULIN AGAR BASE (BBE AGAR BASE). Mix well and pour in sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

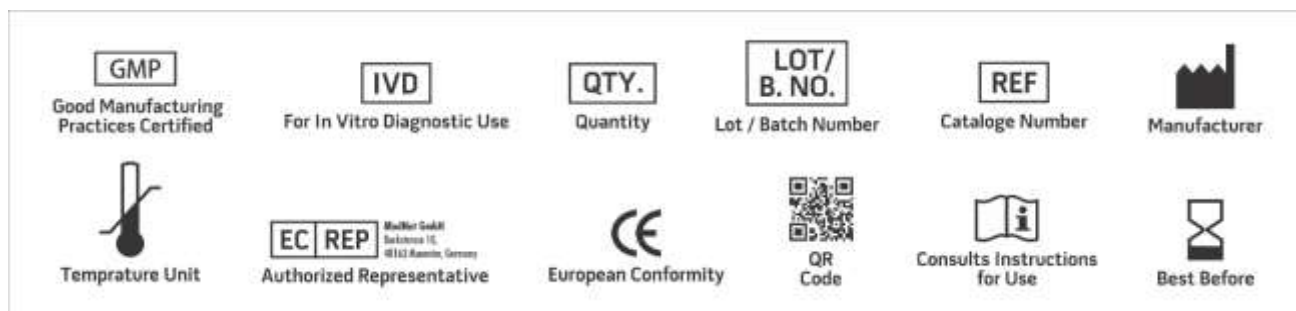

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 5 March., 2022