

TS 116 – LEGIONELLA BMPA SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of *Legionella* species.

COMPOSITION

Ingredients	Concentration
Cefamandole	0.400mg
Polymyxin B	8000IU
Anisomycin	8mg

(Per vial sufficient for 500 ml medium)

Applied in Medium

TM 1218 - LEGIONELLA AGAR BASE

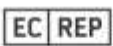
TM 2156 - LEGIONELLA AGAR BASE W/O CHARCOAL

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml of 50% ethanol. Mix well and aseptically add to 450 ml of sterile, molten, cooled (45-50°C) TM 1218 - LEGIONELLA AGAR BASE / TM 2156 - LEGIONELLA AGAR BASE W/O CHARCOAL along with 50 ml of sterile Legionella Growth Supplement (BCYE). Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative MedWor GmbH Buckhorn 10, 48147 Münster, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**

Revision: 7 March., 2022