

TS 128 – GENTA-OXY SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of Vancomycin Resistant Enterococci.

COMPOSITION

Ingredients	Concentration
Oxytetracycline	50mg
Gentamicin	20mg

(per vial sufficient for 500 ml medium)

Applied in Medium

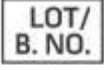
TM 1061 - OXYTETRA GLUCOSE YEAST AGAR BASE (OGYE AGAR BASE) ISO 6611:2004

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml of sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) TM 1061 - OXYTETRA GLUCOSE YEAST AGAR BASE (OGYE AGAR BASE) ISO 6611:2004. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 28 Feb., 2022