

TS 132 – POTASSIUM CHLORATE SUPPLEMENT

INTENDED USE

For selective enrichment of *Yersinia enterocolitica*.

COMPOSITION

Ingredients	Concentration
Potassium chlorate	1g
Distilled water	10ml

(per vial sufficient for 1000 ml medium)

Applied in Medium




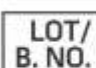







TM 1357 - ITC BROTH BASE (TTC BROTH BASE) (IRGASAN TICARCILLIN AND POTASSIUM CHLORATE BROTH BASE)

INSTRUCTION FOR USE

Warm up the refrigerated solution to room temperature. Aseptically add the contents of 1 vial to 988 ml of sterile, cooled (45-50°C) TM 1357 - ITC BROTH BASE (TTC BROTH BASE) (IRGASAN TICARCILLIN AND POTASSIUM CHLORATE BROTH BASE) along with the contents of 1 vial of Ticarcillin Supplement. Mix well and dispense into sterile tubes.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative <small>MediStar GmbH Buckhorn 16, 48147 Münster, Germany</small>	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 4 March., 2022