

TS 133 – TRICHOMONAS SELECTIVE SUPPLEMENT II

INTENDED USE

For selective isolation of *Trichoderma* species.

COMPOSITION

Ingredients	Concentration
Streptomycin	250mg
Penicillin G	500000Unit

(per vial sufficient for 500 ml medium)

Applied in Medium


TM 1307 - TRICHOMONAS AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml of sterile distilled water. Mix well, avoiding frothing of the solution. Aseptically add the contents of 1 vial to 460 ml of sterile, molten, cooled (45-50°C) TM 1307 - TRICHOMONAS AGAR BASE along with 40 ml of inactivated Horse Serum. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

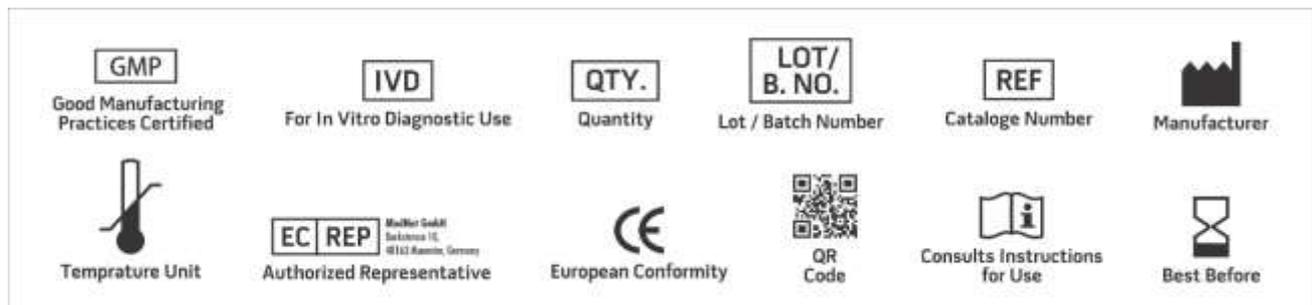

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 3 March., 2022