

TS 152 – MKTT NOVOBIOCIN SUPPLEMENT

INTENDED USE

For selective enrichment and isolation of Salmonellae.

COMPOSITION

Ingredients	Concentration
Novobiocin	40mg

(Per vial sufficient for 450 ml medium)

Applied in Medium




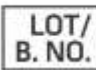








TM 1399 - MULLER KAUFFMAN TETRATHIONATE NOVOBIOCIN BROTH BASE (ISO 6579-1:2017)

INSTRUCTION FOR USE

Rehydrate contents of 1 vial aseptically with 5 ml of sterile distilled water and aseptically add to sterile, cooled (45-50°C) TM 1399 - MULLER KAUFFMAN TETRATHIONATE NOVOBIOCIN BROTH BASE (ISO 6579-1:2017). Mix well and dispense as desired.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 IVD For In Vitro Diagnostic Use	 QTY. Quantity	 LOT/ B. NO. Lot / Batch Number	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 EC REP Authorized Representative <small>Mueller GmbH Südstraße 11 48153 Münster, Germany</small>	 European Conformity	 QR Code	 Consults instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**

Revision: 7 March., 2022