

TS 157 – ENTEROCUCUS SELECTIVE SUPPLEMENT

INTENDED USE

For enumeration of intestinal Enterococci from surface and waste water.

COMPOSITION

Ingredients	Concentration
Nalidixic acid	25mg
Thallos acetate	200mg

(per vial sufficient for 100 ml medium)

Applied in Medium

TM 1388 MUG SF BROTH BASE


TM 2220 MUD SF BROTH BASE

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 5 ml of sterile distilled water. Mix well and aseptically add it to sterile, cooled (45-50°C) 100 ml of TM 1388 - MUG SF BROTH BASE / TM 2220 - MUD SF BROTH BASE along with 1 ml of sterile 1% T.T.C. Solution. Mix well and dispense as desired.

STORAGE

Vials should be stored in sealed container at 2-8°C.

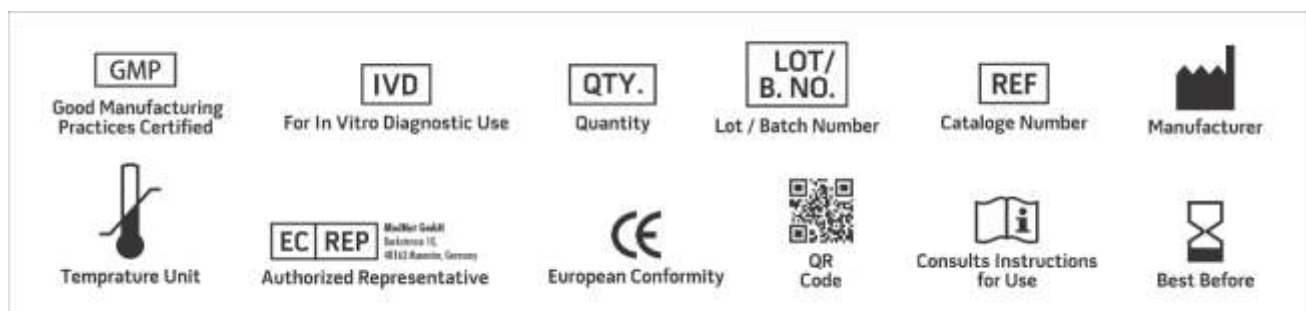

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 28 Feb., 2022