

TS 159 – PERFRINGEN'S SUPPLEMENT-B

INTENDED USE

For selective isolation of *Clostridium perfringens*.

COMPOSITION

Ingredients	Concentration
Polymyxin B sulphate	5000IU
Oleandomycin phosphate	0.250mg

(per vial sufficient for 500 ml medium)

Applied in Medium

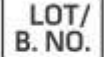
TM 808 - PERFRINGENS AGAR BASE (O.P.S.P)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2 ml sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) TM 808 - PERFRINGENS AGAR BASE (O.P.S.P) along with sterile rehydrated contents of 1 vial of Perfringens Supplement A. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 4 March., 2022