

## TS 184 – GBS SUPPLEMENT

### INTENDED USE

For selective isolation and presumptive identification of group B Streptococci.

### COMPOSITION

Ingredients	Concentration
Metronidazole	5mg
Gentamicin	1mg

(per vial sufficient for 500 ml medium)

### Applied in Medium


TM 1545 - GBS MEDIUM BASE

### INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2 ml of sterile distilled water. Mix well and add it to 475 ml sterile, molten, cooled (45-50°C) TM 1545 – GBS MEDIUM BASE along with 25 ml inactivated Horse serum TS 014. Mix well.

### STORAGE

Vials should be stored in sealed container at 2-8°C.

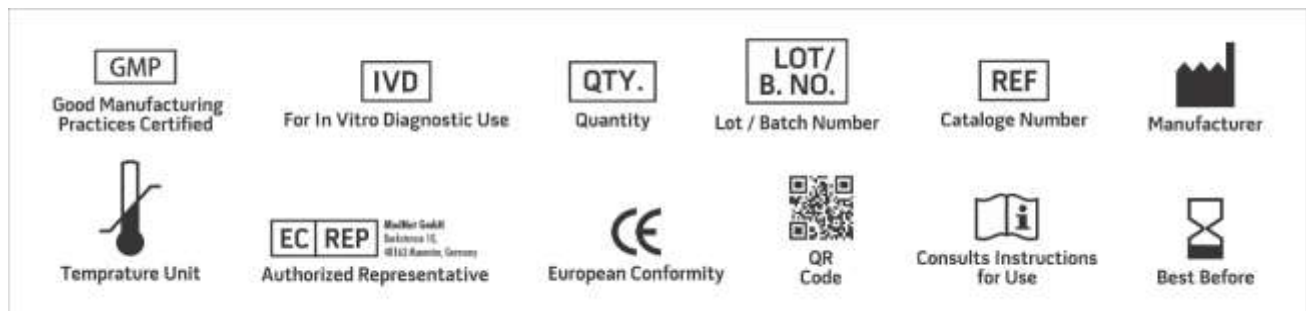

**NOTE:** Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

**\*For Lab Use Only**  
Revision: 28 Feb., 2022