

TS 213 – CHROMOGENIC CANDIDA SELECTIVE SUPPLEMENT

INTENDED USE

For isolation and identification of *Candida* species from mixed cultures.

COMPOSITION

Ingredients	Concentration
Gentamicin	50mg

(Per vial sufficient for 500 ml medium)

Applied in Medium


TM 1588 - CHROMOGENIC CANDIDA DIFFERENTIAL AGAR MODIFIED

INSTRUCTION FOR USE

Rehydrate the content of 1 vial with 5 ml of sterile distilled water. Mix well and aseptically add to 500 ml of sterile, molten, cooled (45-50°C) TM 1588 - CHROMOGENIC CANDIDA DIFFERENTIAL AGAR MODIFIED. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

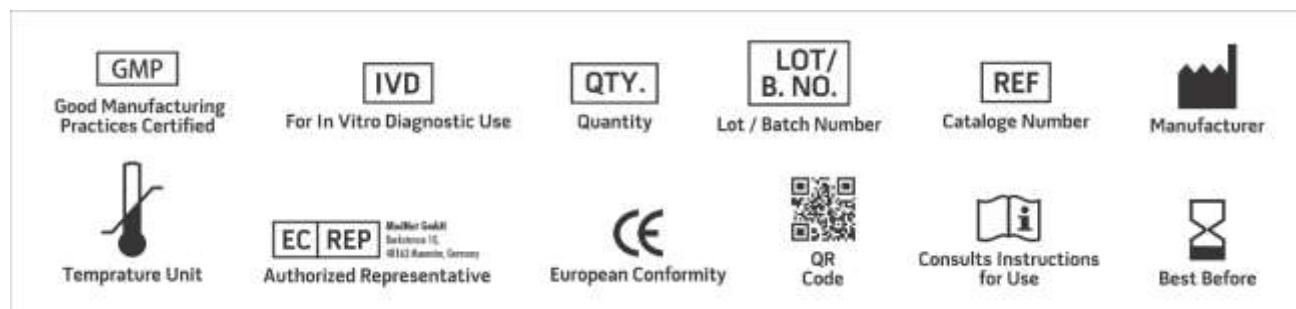

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 5 March., 2022