

TS 229 – L.MONO ENRICHMENT SUPPLEMENT II

INTENDED USE

For selective differentiation of *Listeria monocytogenes* from other *Listeria* species.

COMPOSITION

Ingredients	Concentration
L-phosphatidylinositol	0.5g
Distilled water	15ml

(per vial sufficient for 500 ml medium)

Applied in Medium


TM 1806 - L. MONO CONFIRMATORY AGAR BASE

INSTRUCTION FOR USE

Thaw the contents of 1 vial of L. mono Enrichment Supplement II at room temperature. Aseptically add the sterile contents of one vial to 470 ml of sterile, molten, cooled (45-50°C) TM 1806 - L. MONO CONFIRMATORY AGAR BASE along with sterile rehydrated contents of one vial each of L. mono Selective Supplement I (TS 227) and L. mono Selective Supplement II (TS 228). Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

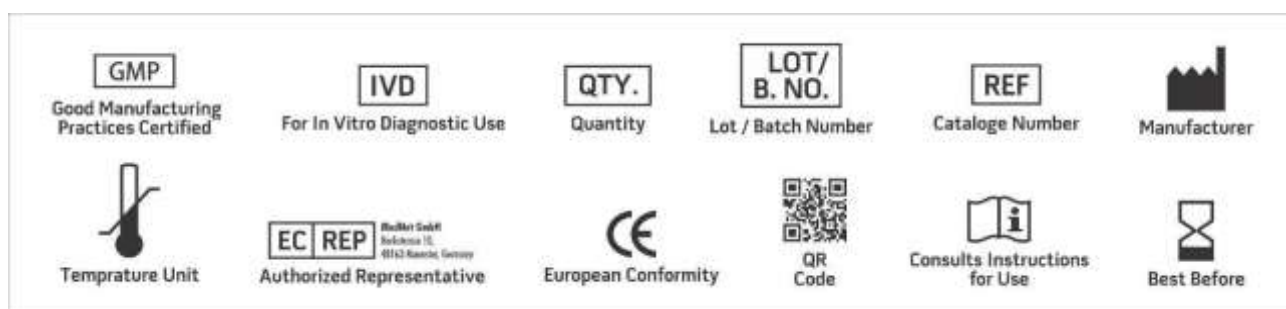

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 05 Aug., 2023