

## TS 231 – MeReSa SELECTIVE SUPPLEMENT

### INTENDED USE

For selective isolation of Methicillin Resistant *Staphylococcus aureus* from clinical specimens.

### COMPOSITION

Ingredients	Concentration
Methicillin	2mg

(Per vial sufficient for 500 ml medium)

### Applied in Medium


TM 1561 - MeReSA AGAR BASE

### INSTRUCTION FOR USE

Rehydrate the contents of one vial with 5 ml of sterile distilled water and aseptically add to 500 ml of sterile, molten, cooled (45-50°C) TM 1561 - MeReSA AGAR BASE. Mix well and pour into sterile petri plates.

### STORAGE

Vials should be stored in sealed container at 2-8°C.

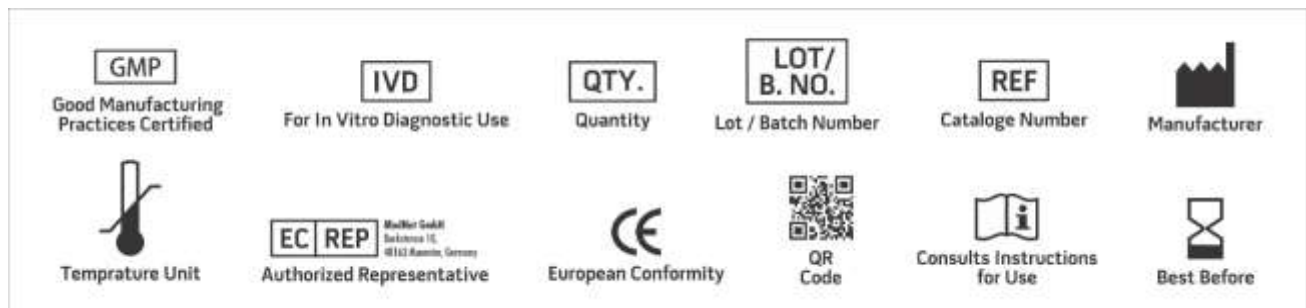

**NOTE:** Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

**\*For Lab Use Only**  
Revision: 7 March., 2022