

TS 252 – SALMONELLA SELECTIVE SUPPLEMENT

INTENDED USE

For isolation and identification of *Salmonella* spp.

COMPOSITION

Ingredients	Concentration
Novobiocin	2.5mg
Cefsulodin	6mg

(per vial sufficient for 500 ml medium)

Applied in Medium


TM 1839 - CHROMOGENIC SALMONELLA AGAR

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5ml of sterile distilled water. Mix gently to dissolve the contents completely. Aseptically add the rehydrated contents of one vial to 500ml of sterile, cooled (45-50°C) of TM 1839 - CHROMOGENIC SALMONELLA AGAR. Mix well and dispense as desired.

STORAGE

Vials should be stored in sealed container at 2-8°C.

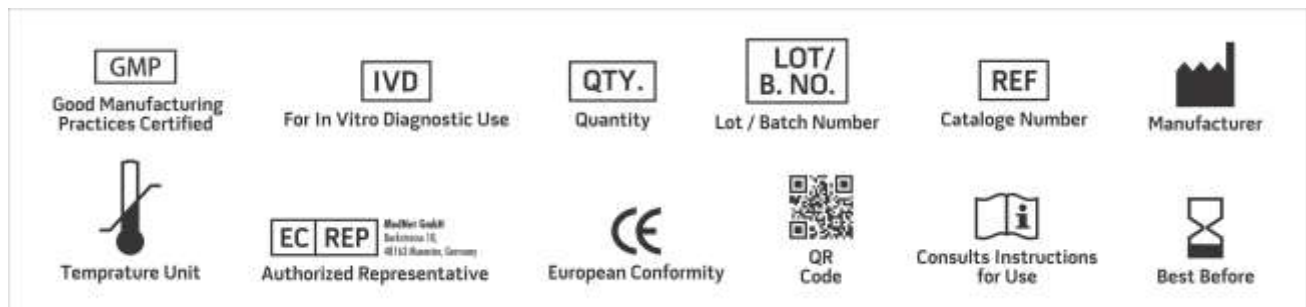

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 3 March., 2022