

## TS 287 – ECC SELECTIVE SUPPLEMENT MODIFIED

### INTENDED USE

Recommended for detection of *Escherichia coli* and other coliform.

### COMPOSITION

Ingredients	Concentration
Cefsulodin	10mg

(Per vial sufficient for 500 ml medium)

### Applied in Medium


TM 2127 - CHROMOGENIC M-COLICONFIRM BROTH

### INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 2 ml sterile distilled water. Mix well and aseptically add to 1000 ml of sterile, molten, cooled (45-50°C) TM 2127 - CHROMOGENIC M-COLICONFIRM BROTH. Mix well and pour into sterile petri plates.

### STORAGE

Vials should be stored in sealed container at 2-8°C.

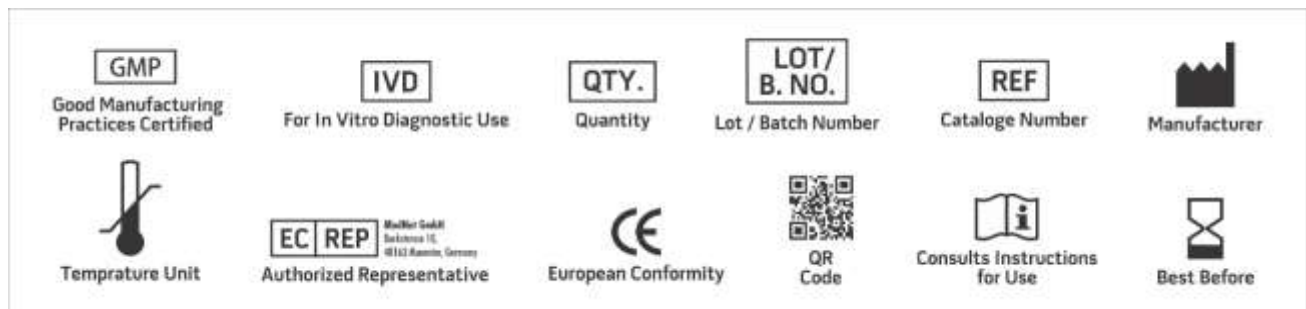

**NOTE:** Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

**\*For Lab Use Only**  
Revision: 11 March., 2022