

TS 318 – M-CP SELECTIVE SUPPLEMENT I

INTENDED USE

Recommended for selective isolation of *Clostridium perfringens*.

COMPOSITION

Ingredients	Concentration
D-Cycloserine	200mg
Polymyxin B sulphate	12.500mg

(Per vial sufficient for 500 ml medium)

Applied in Medium


TM 2130 - M-CP AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 5 ml of sterile distilled water. Mix well and aseptically add it to 485 ml sterile, molten, cooled (45-50°C) TM 2130 - M-CP AGAR BASE along with one vial of M-CP Selective Supplement II / M-CP Selective Supplement, Modified. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

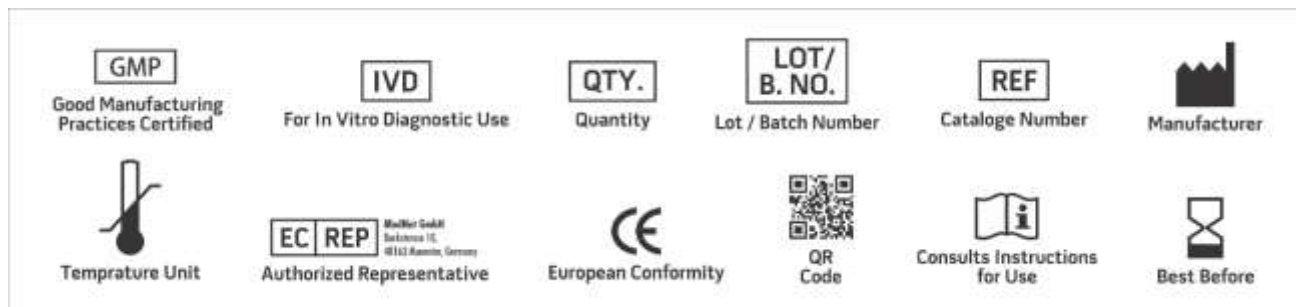

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 7 March., 2022